


IMPORTANT: About Your New AirSense™ 10 Machine

It is important you understand this new device has the ability to transfer data to the manufacturer, ResMed. The airplane mode has been set to "On" to prevent data from being transmitted. When you turn on your machine, a notice will appear like the one shown to the right. This message will ask you to confirm your settings.


If you select "No" your data **will not** transfer to ResMed.

If you select "Yes" then your data **will** transfer to ResMed.


Please Note: Due to HIPAA regulations, your patient information (name, contact information, etc) has not been registered with ResMed.

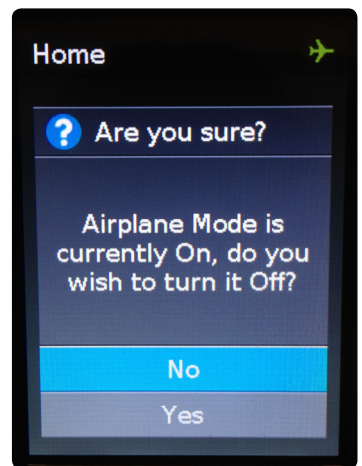


Turning Airplane Mode On for the AirSense™ 10 Line to Prevent Data Transmission

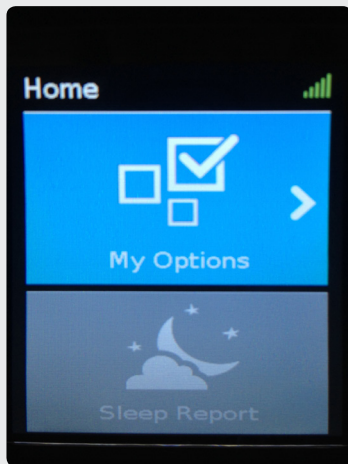


1. On the User Menu, click on "My Options"

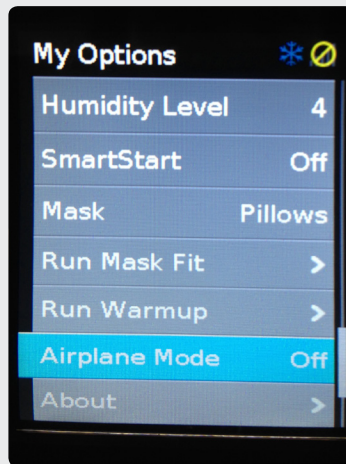

2. Towards the bottom select "Airplane Mode"

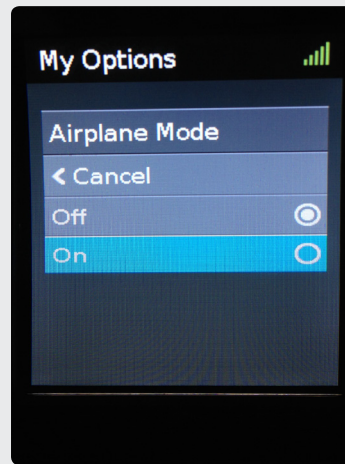

3. Change the default from "Off" to "On"

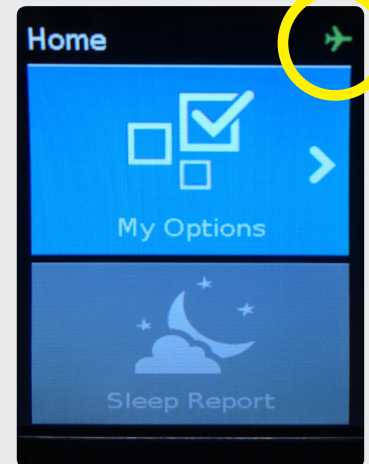

4. Once complete, the airplane icon will show


Machine Icons and Data Transfer



The bar signal icon on the top right hand side of the screen indicates airplane mode is "Off" and the modem is turned on. This setting enables cellular signal and transfers data to the manufacturer of the machine, ResMed.



The circle with line through icon on the top right hand side of the screen indicates the modem is turned on but that there is no cellular signal currently being received. This means airplane mode is "Off" and data can be sent to ResMed once signal is established.

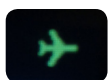

The airplane icon on the top right hand side of the screen indicates airplane mode is "On" and the modem is turned off. This setting prevents data from being transferred to ResMed.

